

Recommended Expansion Protocol for RoosterVial-hUC-1M-XF

Protocol Summary


To expand one vial of xeno-free, human umbilical cord-derived Mesenchymal Stem/Stromal Cells (RoosterVial-hUC-1M-XF) to at least 20 million cells within one week you will need the following reagents, cell culture materials, and equipment.

Materials & Equipment

| Item | Quantity | Vendor | Part Number* |
|---|--------------|-------------------------------|--------------|
| RoosterVial-hUC-1M-XF | 1 Vial | RoosterBio | C43001UC |
| RoosterNourish™-MSC-XF | 1 Bottles | RoosterBio | K82016 |
| Either of the following: | | | |
| T75 CellBIND flasks | 4 | Corning | 3290 |
| T225 CellBIND flasks | 2 | | 3293 |
| TrypLE Select Enzyme | 20 mL | Life Technologies | 12563029 |
| Vitronectin Recombinant Human Protein (500 µg/ml) | ThermoFisher | 1 mL: A14700 10 mL: A31804 | |
| DPBS without Calcium, Magnesium (-/-) | ThermoFisher | 14190144 | |
| Biosafety Cabinet | | | |
| Centrifuge | | | |
| Incubator | | | |
| Water Bath (or ThawSTAR®) | | | |

Process Overview

Cell Expansion Summary

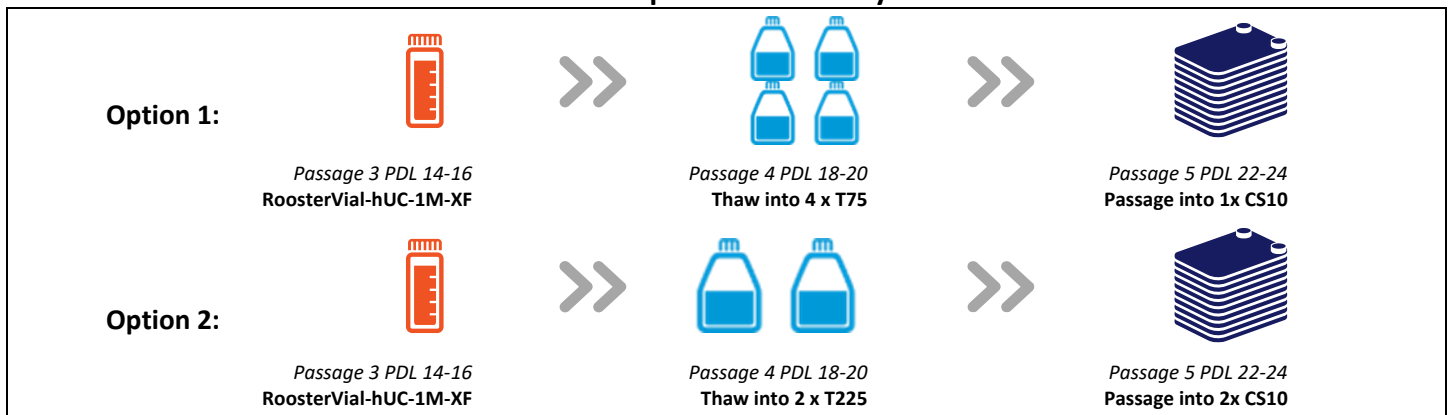

**RoosterBio strongly recommends the use of CellBIND with Vitronectin or Fibronectin surfaces for expansion of the Xeno-Free hUC product line. **Please refer to protocol for full process instructions.*

CONFIDENTIAL. FOR SOLE USE OF RECIPIENT. ALL PRODUCTS ARE FOR RESEARCH USE ONLY.

- Thaw and seed cells at recommended: 3,000 cells/cm² (min. >2,000 cells/cm²).
- Expand cell cultures 3-6 days to >80% confluency at 37°C, 5% CO₂ incubation.
- **NO MEDIA EXCHANGES REQUIRED.** RoosterNourish-MSX-XF does not need to be exchanged, or fed, within 6 days of flask-based culture.

Recommended Protocol

1. Expansion Options

Please refer to the following Expansion Options table to determine the cell culture vessel best suited to your research needs. Note: If larger cell numbers are required, a decrease in seeding density to 2,200 cells/cm² (and increase in total vessel seeding surface area to 450 cm² and total seeding media to 90 mL) will maximize total cell output at harvest with an increase in total Population Doubling Level (PDL).

| Vessel | Surface Area (cm ²) | Number Vessels Needed | Total Surface Area (cm ²) | Seeding Density (cells/cm ²) | Approx. Yield at Harvest | Days of Culture |
|--------|---------------------------------|-----------------------|---------------------------------------|--|--------------------------|-----------------|
| T75 | 75 | 4 | 300 | 3,300 | ≥ 20M | 3 to 6 |
| T225 | 225 | 2 | 450 | 2,200 | ≥ 30M | 3 to 6 |

2. Media Preparation

- 2.1. Bring RoosterNourish-MSX-XF components to room temperature, protected from light, for up to four hours.
- 2.2. Prepare 1 bottle of medium by aseptically adding 1 bottle of RoosterBooster™-MSX-XF (Part No. SU-016) to 1 bottle of RoosterBasal™2.0-CC (Part No. M22520).
- 2.3. Mix well by capping and gently mixing the bottle.

3. Flask Preparation

- 3.1. Upon receipt, thaw the vial of vitronectin and transfer DPBS -/- at room temperature.
- 3.2. Calculate the working concentration of vitronectin (target 0.15 µg /cm²):

| Culture Surface Area (=A) | | | | | | | | |
|---------------------------------|---------|---------|--------|-----|-----|------|-------------|--------------|
| | 24-Well | 12-Well | 6-Well | T25 | T75 | T225 | CellStack 5 | CellStack 10 |
| Surface Area (cm ²) | 2 | 4 | 10 | 25 | 75 | 225 | 3180 | 6360 |

| Raw Data | | Adjusted Data | |
|-------------------|---------------------------|------------------------------|---|
| Surface Area (=A) | Target Concentration (=B) | Total µg VN Required (C)=A*B | # mL of VN Stock Solution (D)=C/500 µg/ml |
| | 0.15 µg/cm ² | | |

| Surface Area (=A) | Volume Conversion (=E) | Total Volume of DPBS Required (F)=A*E |
|-------------------|------------------------|---------------------------------------|
| | 0.1 mL/cm ² | |

- 3.3. Prepare working aliquot of Vitronectin in DPBS -/-.
 - 3.3.1. Aliquot required volume of Vitronectin (**D**).
 - 3.3.2. Dilute in appropriate volume of DPBS -/- (**F**).
 - 3.3.3. Gently resuspend by pipetting the vitronectin dilution up and down.

CONFIDENTIAL. FOR SOLE USE OF RECIPIENT. ALL PRODUCTS ARE FOR RESEARCH USE ONLY.

- 3.3.4. Aliquot remaining Vitronectin and return to -80°C promptly.
- 3.4. Add total volume (**F**) to each vessel.
- 3.5. Incubate flasks at room temperature for at least 1 hour (away from direct light).

| |
|-----------------------|
| Total Incubation time |
| |

- 3.6. Aspirate volume within culture flasks (do not allow vessels to completely dry).

4. Cell Thawing & Seeding

- 4.1. Aseptically transfer 10 mL of prepared medium into a 50 mL centrifuge tube.
- 4.2. Thaw RoosterVial-hUC-1M-XF vial in an automated thawing device (e.g., ThawStar), or manually in a 37°C water bath. When thawing in a water bath, monitor the vial closely and remove from water bath once only a small bit of ice is remaining (2-3 min).
- 4.3. Aseptically transfer vial into a Biosafety Cabinet (BSC).
- 4.4. Transfer vial contents into the 50 mL centrifuge tube containing prepared medium and mix cell suspension well.
- 4.5. Centrifuge at 350 x g for 6 min.
- 4.6. Aspirate the supernatant and resuspend cells in 20 mL of RoosterNourish-MSX-XF medium.
- 4.7. Mix well and seed cells equally into four T75 vessels or into two T225 vessels, and add medium to bring volume up to final volume according to table below:

| Type of Culture Vessel | Total Volume of Cell Suspension per Vessel | Final Volume per Vessel |
|------------------------|--|-------------------------|
| T75 x 4 | 5 mL | 15 mL |
| T225 x 2 | 10 mL | 45 mL |

- 4.8. Transfer vessels into an incubator (37°C, 5% CO₂) and ensure surfaces are covered evenly with media.

5. Cell Expansion

- 5.1. Microscopically monitor cell confluency starting on day 3 of culture.
- 5.2. When culture is >80% confluent, cells are ready to harvest.

| Day | 3 | 4 | 5 | 6 |
|------------------------|---|---|---|---|
| Cell Confluency | | | | |

Note: For best expansion and functional performance, it is recommended to passage the cultures before reaching 90% confluence. If the cultures reach over confluence, this may result in difficulty when harvesting, increased aggregation, decreased cell viability, growth inhibition and loss of differentiation potential.

6. Cell Harvest

- 6.1. For harvest, transfer vessels into biosafety cabinet and remove spent media.
- 6.2. Add 3 mL TrypLE to each T75 flask or 10 mL TrypLE to each T225 flask.
- 6.3. Distribute TrypLE evenly to cover all the cells and place vessels in 37°C (5% CO₂) incubator. Check culture every 5 min until cells are detached from surface. Gently tap to dislodge remaining cells from surface.

CONFIDENTIAL. FOR SOLE USE OF RECIPIENT. ALL PRODUCTS ARE FOR RESEARCH USE ONLY.

Total Time Required for Cell Detachment

- 6.4. Add equivalent volume of RoosterNourish™-MSC-XF to each vessel to stop the TrypLE activity.
- 6.5. Transfer the cell suspension into a 50 mL centrifuge tube.
- 6.6. Centrifuge at 350 x g for 6 min.
- 6.7. Aspirate the supernatant.
- 6.8. Resuspend cells in medium to achieve desired cell concentration. Measure the total volume of cell suspension:

Total Volume of Cell Suspension (=A)

- 6.9. Transfer 0.5 mL of cells into microcentrifuge tubes for cell counts.
- 6.10. Count cells with a cell counting device, performing a dilution if required to get within its acceptable range:

| Raw Data | | Adjusted Data | |
|----------------------|--------------------------------|----------------------------|--------------------------------|
| Dilution Factor (=B) | Viable Cell Concentration (=C) | Cell Concentration (D)=B*C | Total Cells at Harvest (E)=D*A |
| | | | |

- 6.11. Cells are ready to be used in your application.

7. Expansion Options —2nd Passage

Note: One passage expansion is described above. If additional cells are required, second passage expansion options are listed below (select option based on required cell number at harvest). It is recommended that cells be used within 10-13 PDLs of the working cell bank.

| Vessel | Surface Area (cm ²) | Number Vessels Needed | Total Surface Area (cm ²) | Seeding Density (cells/cm ²) | Approx. Yield at Harvest | Days of Culture |
|--------|---------------------------------|-----------------------|---------------------------------------|--|--------------------------|-----------------|
| CS10 | 6360 | 1 | 6,360 | 3,000-6,000 | ≥ 500M | 3 to 6 |
| CS10 | 6360 | 2 | 12,720 | 2,200-4,200 | ≥ 700M | 3 to 6 |

Caution to Users: RoosterBio products contain human sourced materials and should be treated as potentially infectious. Employ universal safety precautions and wear protective clothing and eyewear while handling. Practice appropriate disposal techniques per CDC guidelines for biohazardous material.

Provision of Seller Product subject to Seller Standard Terms and Conditions. Any technical advice furnished, or recommendation made concerning any use or application of any Seller Product is believed to be reliable, but Seller makes no warranty, either express or implied, as to its accuracy or completeness or of the results to be obtained. Purchaser assumes full responsibility for quality control, testing and determination of suitability of Seller Product for its intended application or use.

CONFIDENTIAL. FOR SOLE USE OF RECIPIENT. ALL PRODUCTS ARE FOR RESEARCH USE ONLY.